


Installations Manual
ParisServer

1. Contents

1. Contents.....	2
2. Installation.....	3
3. Configuration SQL.....	4
3.1 SQLManager.....	4
4. Configuration VarioStorage.....	5

2. Installation

- Click “Next”

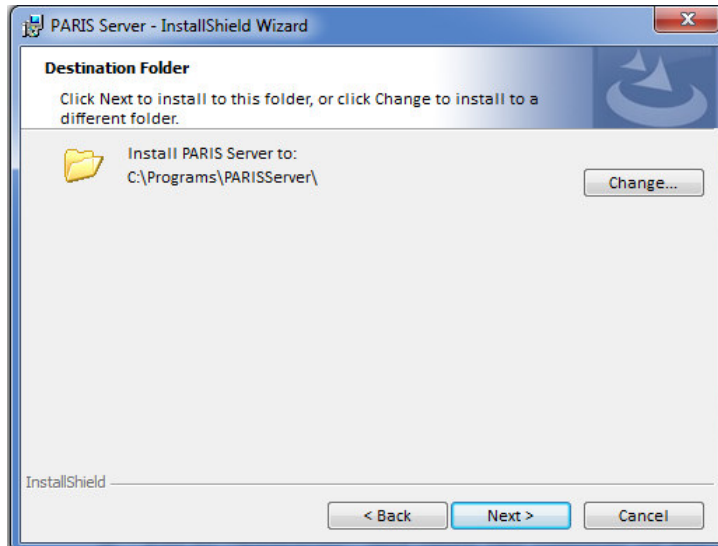


- Choose as Path e.g. “C:\Programs\PARISServer”
- Click “Next”
- Click “Install”


3. Configuration SQL

3.1 SQLManager

- Open SQLManager

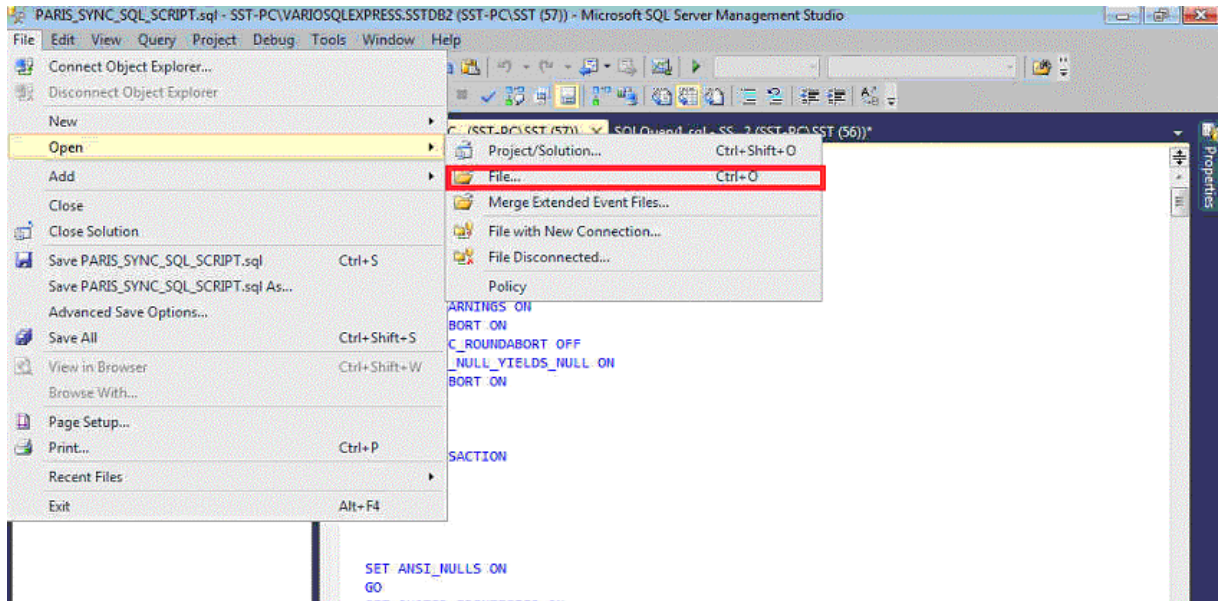


- Click on “File”, go to “Open” → “File” and search for the File “PARIS_SYNC_SQL_SCRIPT.sql” inside the installation folder of PARISServer. (E.g. C:\Programs\PARISServer\SQL\PARIS_SYNC_SQL_SCRIPT.sql)

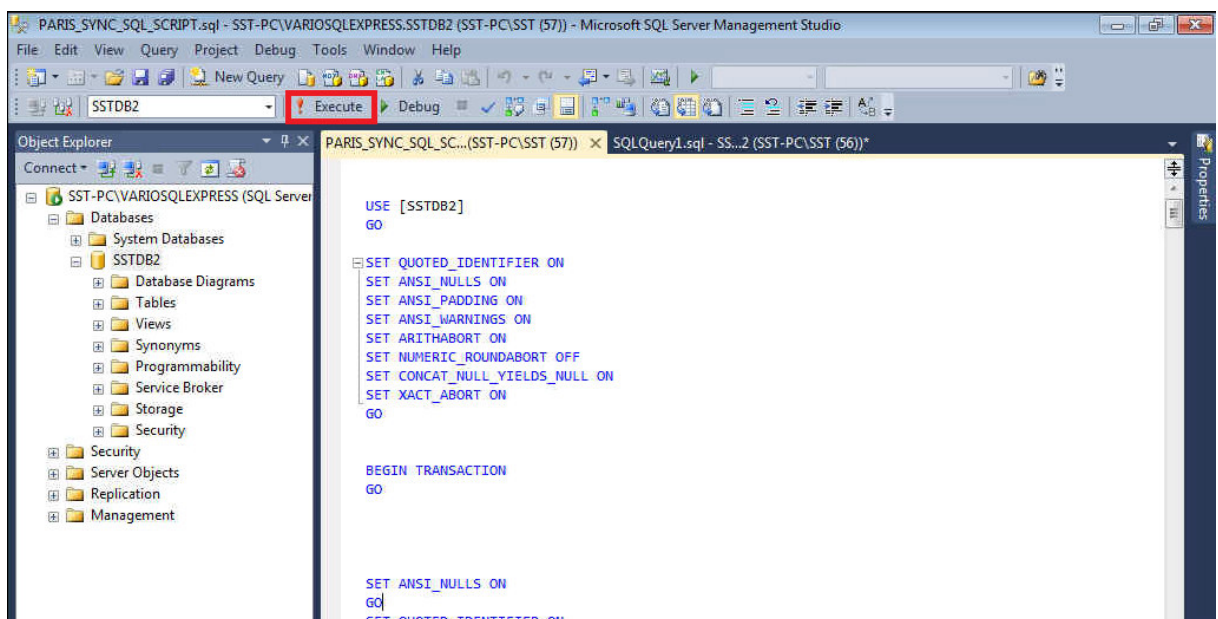

- After opening the file click “Execute”

4. Configuration VarioStorage

- Open VarioStorage
- Go to “Services” and open the settings of ParisServer

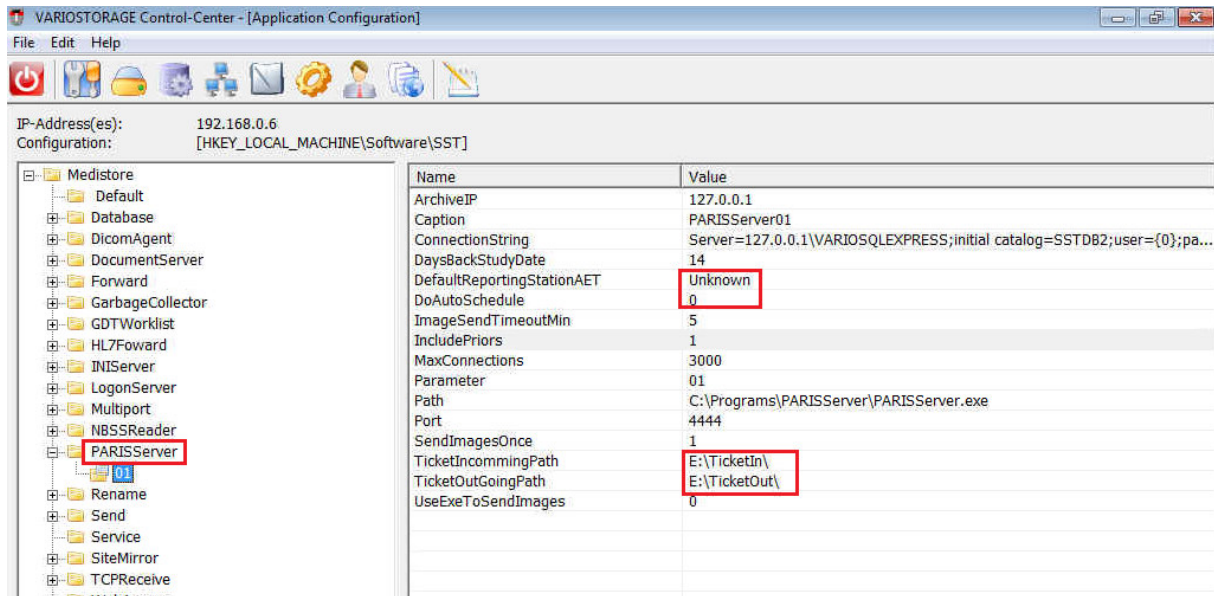

- Insert in “DefaultReportingStationAET” : “Unkown”
 - Insert in „TicketIncommingPath“ : „E:\TicketIn“
 - Insert in “TicketOutGoingPath” : “E:\TicketOut”
 - Insert in “DoAutoSchedule” : “0”
-
- Start the Service